

## Data sheet


## Fume hood Compact Line DCL-1200 CR



## Standard equipment for fume hoods:


- structure entirely made of galvanized steel without the use of wood-based materials
- worktop at a height of 900mm made of monolithic ceramics with a raised rim with a suspended ceramic sink 280x80mm
- toughened glass window, lifted manually with counterweight, opened upwards with a lock at the height of 500 mm (max 810mm), and additionally sliding glass (right / left). The window has a system preventing the uncontrolled drop of the window
- 2 electric sockets 230V 50Hz
- 2 water spouts with valves on the BROEN front panel
- energy-saving linear LED lighting separated from the fume hood with a glass pane
- double system the back wall with exhaust slots in the upper and lower part ensuring proper extraction of vapors lighter or heavier than air without turbulence and dead zones
- safety flap in the ceiling ensuring decompression in the event of uncontrolled pressure increase in the working chamber
- full monitoring set Schneider brand) air flow velocity in the fume cupboard window compliant with the PE-EN 14175 standard
- CE Declaration of Conformity: PN-EN 14175; PN-EN 14056; PN-EN 13150



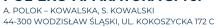





TECHNICAL DATA - DIMENSIONS		
overall dims [mm]		
width A	1280	
height B	2325-2600	
depth C	960	
working space dims (mm)		
width D	1150	
window width D1	965	
height E	1220	
depth F	635	

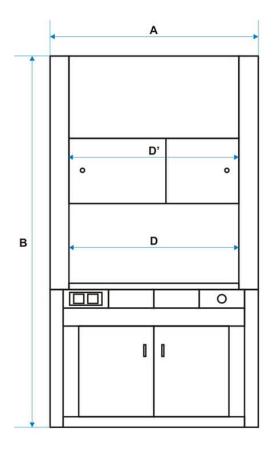

tel. (+48) 32 453 91 70

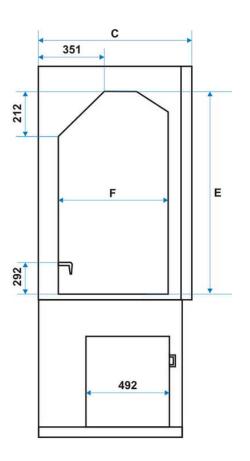

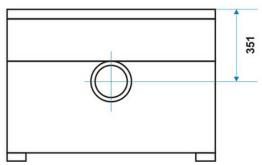



A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

ecommended airflow [m3/h]	600 - 950
equired capacity of the ventilation system [m/s]	0,3 - 0,5
ominal power [W]	46
roltage**	230 VAC 50Hz
liameter of the air outlet stub [mm]	Ø 200
vater supply	G 1/2 " male thread
ewage connection diameter [mm]	Ø 50
electrical insulation	class I
tandard ventilation system	double rear wall
llumination installed outside the working chamber	LED, class A ++, through insulating glass / from operator panel
nonitoring set	standard: - FM 550-A-0-E controller with a speed sensor AFS 100 in the fume hood window and the FA-0025-3 operator panel option: - ICM 500-F-EU controller with flow control - ICM 500-FP-EU controller with flow control and cooperation with BMS - other Schneider controllers available on request
aising the window	standard: - manual with counterweight (lock at a height of 500mm, maximum opening 810mm) option: - SC 500 - Auto Protect system - with foot switch - SC 500 - Auto Protect system - with panel switch
rame and housing material	The frame is made entirely of galvanized steel. External housing made of galvanized and painted with epoxy paint steel
vorking chamber wall finishing	Buchtal large-size ceramic
vorktop	Monolithic ceramics with marine edge
vorktop OPTION (available materials)	- Phenolic resin, epoxy - Stainless steel AISI 304, 316
inderbench cabinet OPTION (different types)	<ul> <li>metal, ventilated cabinet painted with epoxy paint</li> <li>cabinet for acids and alkalis made of foamed PVC</li> <li>explosion-proof cabinets Asecos 30 and 90 minutes</li> </ul>
varranty	24 months

## Other optional equipment:

- main switch
- anti-shock switch
- safety switch
- 400VAC IP44 socket
- additional 230VAC IP44 sockets
- 230VAC sockets inside the working chamber IP44
- Manual Protect system battery backup of the monitoring set
- power supply control of the socket with the option of timer setting
- additional gas taps
- metal cabinet key lock under the tabletop
- explosion-proof equipment EX
- side glazing
- truss made of acid-resistant steel, on the back wall of the manipulation chamber
- inspection openings
- flexible strip brush

